



FURIX

FVA-20

FURIX ANAESTHESIA VENTILATOR MACHINE



PET HEALTH CHECK

- ✓ Deliver ventilation
- ✓ Comprehensive consumables
- ✓ Support various scenarios
- ✓ Support anaesthesia
- ✓ Comprehensive solution

VENTILATOR SPECIFICATION

Drive	Pneumatically driven and electronically control
Working mode	Standby/Manual/ACGO Volume Support (VS) Volume Control Ventilation (VCV) Pressure Control Ventilation (PCV)
VT	0mL-3000mL
MV	0L/min~100L/min
PEAK	-20cmH2O~120cmH2O
PEEP	0cmH2O~70cmH2O
RR	0bpm-120bpm
VT	<75mL: ±15mL ≥ 75mL: ± 20mL or ± 10% of the set value, whichever is greater
MV	±0.2L/min or ± 15% of the set value, whichever is greater
Ppeak	(2% of full scale + 4% of the readings)
PEEP	OFF: Do not define 4cmH2O~30cmH2O: ±2.0cmH2O or + 10% of the set value, whichever is greater
RR	±1bpm or ± 5% of the set value, whichever is greater
Active AGSS Size	500mm x150mm x120mm
Active AGSS Pump rate	80L/min
Passive AGSS Connector	30mm OD conical 300mL
Small bellow capacity	300mL
Small bellow capacity	1500mL



Anesthetic Vaporizer Specifications

Vaporizer	Pour Fill (Isoflurane / Sevoflurane / Enflurane)
Filling methods	
Weight	8.0kg (empty) 8.5kg (full)
Filling volume	360ml (dry wick) 300ml (moist wick) 260ml (between the minimum and maximum marks)
Concentration range	0-6% (Isoflurane, Enflurane) 0-8% (Sevoflurane)
Concentration accuracy range	-20%~+30% of the set value or -5% ~ +7.5% of the max value, whichever is the greater
Vaporizer base	Single position, Selectatec connecting type

Application Scenarios




Touch-screen, easy to use

10.1" screen color touch screen Intuitive partition interface display



Reforming designs

Multiple mounting solutions for Patient monitors, pumps and oxygen generators

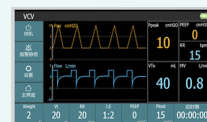

Powerful Ventilation

Multiple ventilation modes Tidal volume: 5ml-1500ml Smart parameters setting



Furix Trade Venture Limited

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,
Mob:+852 9579 0472 E-mail: vet.furix@gmail.com, sales@myfurix.com
Tel:+852 9579 0472 www.myfurix.com

FURIX

Veterinary Anesthetic Ventilator Machine



Touch-screen, easy to use



Reforming designs




Powerful Ventilation



Deliver Ventilation

- Customized alarm upper and lower limits Animal weight, lung condition
- Three-level sound and light alarm Sound frequency, background color prompt

待机	报警限	当前报警
报警静音	下限	上限
设置	Paw cmH2O	4 30
主界面	RR bpm	OFF OFF
	MV L/min	2 6
	VT mL	1 1000
		1 2 3 <x>
		4 5 6 OFF
		7 8 9
		0 . ✓

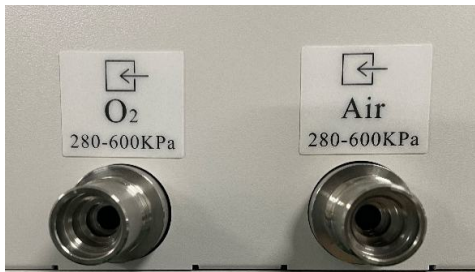

- Air source O2 Air
- Ventilator VCV PCV VS Mini tidal volume 5ml
- Gas ETCO2 O2

Reduce pressure within respiratory management to avoid pulmonary pressure injury
For adjusting the pressure limits of breathing circuit systems

Mechanical Ventilation ✗ Manual Ventilation ✓



Flush O2



Waste canister AGSS



ACGO Switch



ISO Vaporizer

Furix Trade Venture Limited

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,
Mob: +852 9579 0472 E-mail: vet.furix@gmail.com, sales@myfurix.com
Tel: +852 9579 0472 www.myfurix.com

FURIX